





PAL SPME Arrow The Better SPME

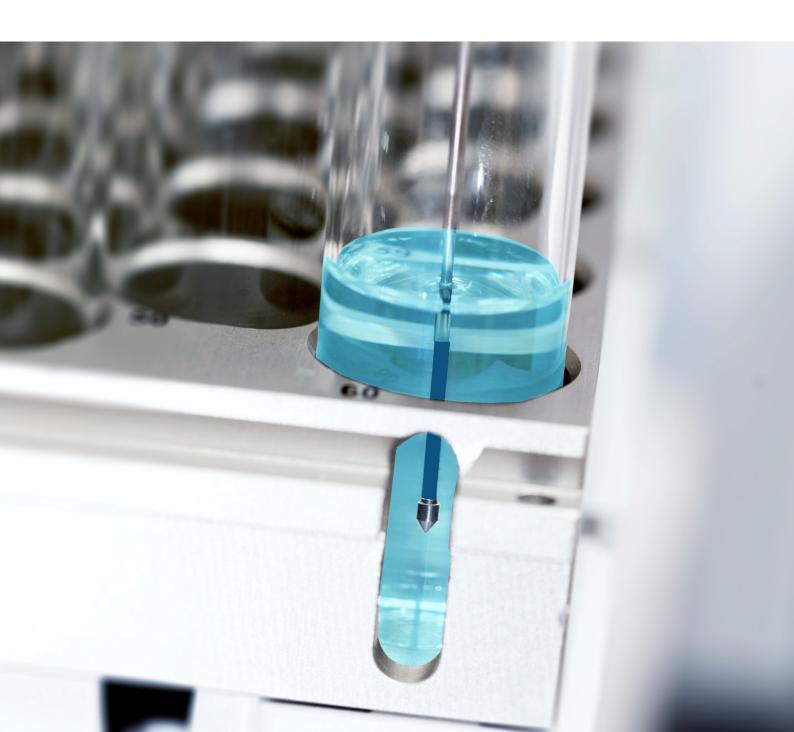

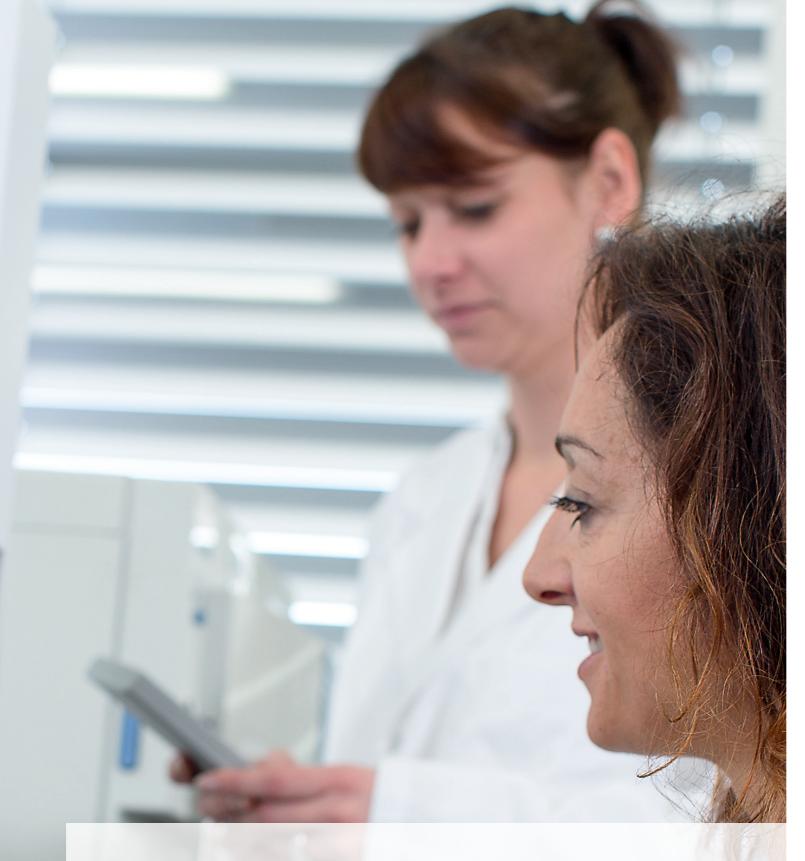


PAL SPME Arrow: The new dimension of SPME

PAL SPME Arrow The new dimension for Solid-Phase Micro Extraction

SPME has become one of the most widely used extraction technologies for environmental, food and clinical analyses. It is well suited for automated sample preparation resulting in reduced time per sample, less sample manipulation and solvent consumption. However, the technology remained almost unchanged with some significant drawbacks, such as the limited mechanical stability and small phase volumes of the fibers.

The PAL SPME Arrow (patent pending) is a new technology for micro-extraction, combining trace level sensitivity with high mechanical robustness. The PAL SPME Arrow has an outer diameter of 1.1 or 1.5mm, resulting in large sorption phase surfaces and volumes. The arrow-shaped tip allows smooth penetration of vial and injector septa. In contrast to traditional SPME fibers, the Arrow design fully protects the sorptive material, minimizing adverse influences and loss of analytes during transfer processes. With the PAL RTC and RSI the SPME Arrow sampling is fully automated leading to high productivity.

The figure below shows the dimension of a PAL SPME Arrow 1.5 mm (a), a PAL SPME Arrow 1.1 mm (b) and a SPME Fiber (c) in comparison.

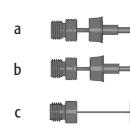




What Better SPME means

Bigger surface, faster extraction. More sorption phase, superior sensitivity. Optimized geometry, greater robustness.

A

M

= 2 x productivity

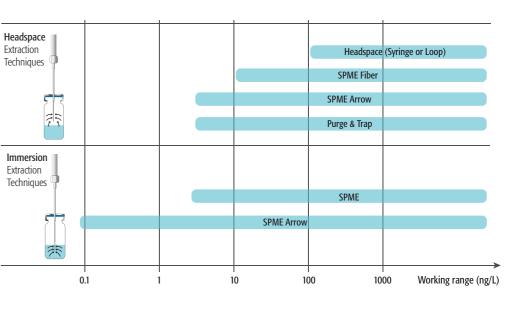

Better SPME

- pounds are now amenable to SPME

PAL SPME Arrow covers a wide range of applications. However, for dynamic headspace applications, especially for volatiles we recommend ITEX Dynamic Headspace (DHS). This powerful technology achieves ng/L sensitivities without the pitfalls of purge & trap systems.

2 x higher sample throughput.

Up to 10 x more sensitivity - wider linear range.

PAL SPME Arrows last at least 2 x longer. Injector septa last at least 2 x longer. Lower running costs.

- Adaptation of existing SPME methods is straightforward


- PAL SPME Arrow works well for headspace and immersion extraction

- With the wide selection of sorption materials (cf. page 10) a wide variety of com-


- The PAL SPME Arrow is an ideal field sampling device

More Volume: Up to 10x more Sensitivity

The table below shows the dimension of a PAL SPME Arrow 1.5 mm (a), 1.1 mm (b) and a SPME Fiber (c) in comparison:

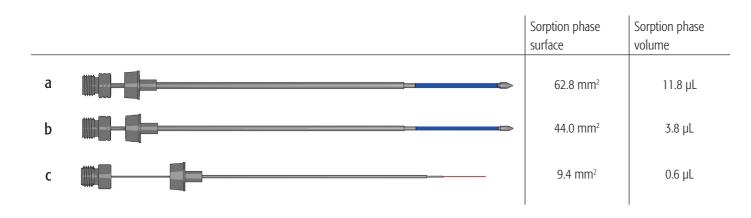


Immersion Extraction: Polyaromatic Hydrocarbons (PAHs) in Water

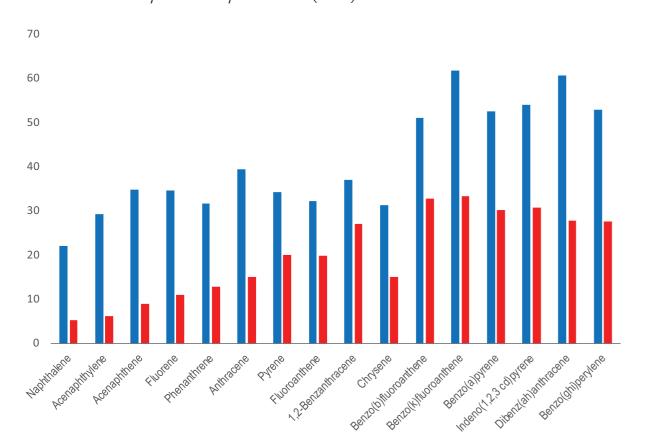


Extraction yields for water samples containing PAHs @ 50 ng/L after 70 min. A PAL SPME Arrow 20 mm x 250 µm Carbon WR was compared to SPME fiber 10 mm x 100 µm Carbon WR immersion (Kremser et al., 2015).

Headspace Extraction: Aroma Analysis in White Wines

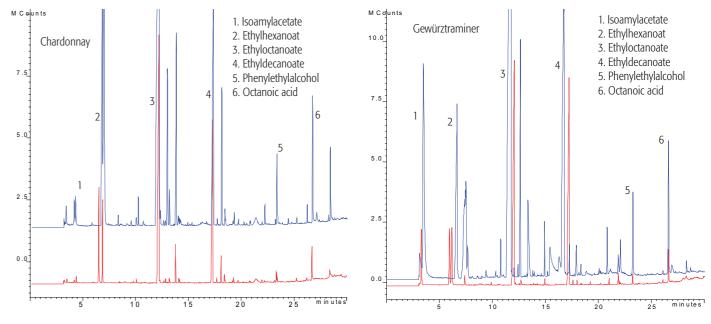


Chromatograms showing the headspace extraction of aroma components from different white wines with PDMS fibers (PAL SPME Arrow 100 µm, 20 x 1.1 mm compared to SPME 100 µm, 10 x 0.3 mm)

Iodoform in Water

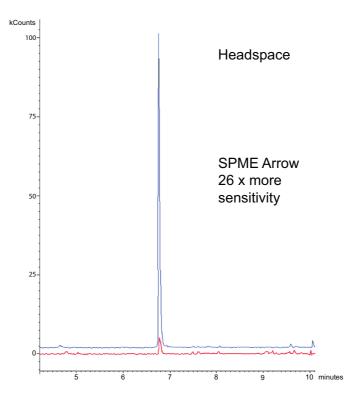


Extraction of 1 µg/L iodoform from tap water with DVB fibers (headspace and immersion extraction), PAL SPME Arrow 100 µm, 20 x 1.1 mm compared to SPME 100 µm, 10 x 0.3 mm.

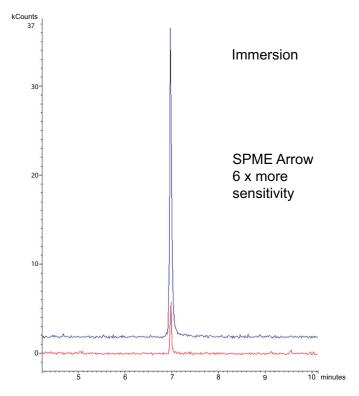


Bigger Surface: 2x Throughput

Immersion Extraction: Polyaromatic Hydrocarbons (PAHs) in Water

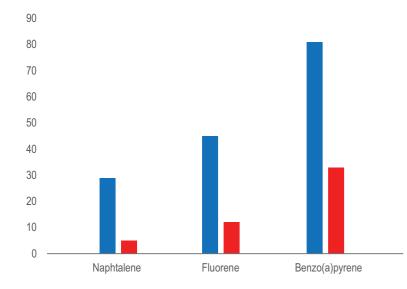

Relative immersion extraction yield (measured as % extracted after a 70 min) for PAHs at 50 ng/L with PDMS fibers (PAL SPME Arrow 100 µm, 20 x 1.1 mm compared to SPME 100 µm, 10 x 0.3 mm)

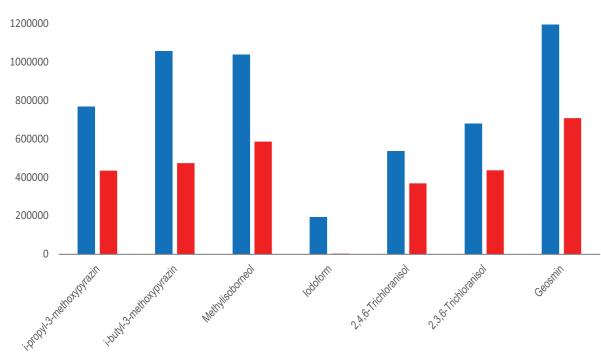

Headspace Extraction: Off Flavor Compounds in Water


DVB SPME Arrow DVB SPME

Relative headspace extraction yield (measured as amount extracted after 30 min) for off-flavor compounds in water at 100 ng/L with DVB fibers. (PAL SPME Arrow 100 µm, 20 x 1.1 mm compared to SPME 100 µm, 10 x 0.3 mm)

References

- [1] Belardi R., Pawliszyn J., Water Pollut. Res.J.Can. 1989, 24, 179
- [2] SPME Arrow Evaluation of a Novel Solid-Phase Microextraction Device for Freely Dissolved PAHs in Water; Kremser A. et al., Anal. Bioanal. Chem. 2016, 408, 943-952
- [3] Solid phase microextraction Arrow for the sampling of volatile amines in wastewater and atmosphere; Helin A. Et al., J. Chrom. A 2015, in press
- [4] PAL System Application Notes: Determination of iodoform in drinking water by SPME and GC/MS and Determination of C2-C12 aldehydes by SPME on-fiber derivatization and GCMS

With the PAL RTC and RSI the entire SPME process is fully automated. This guarantees process safety and high reproducibility.

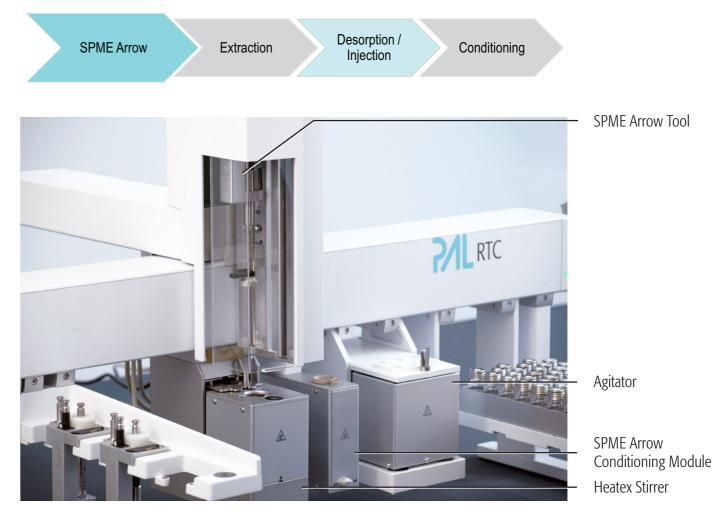

PAL Heatex Stirrer - New Mixing and Heating Technology for Sample Preparation and SPME.

The powerful PAL Heatex Stirrer mixes samples rapidly applying cycloid shaped mixing patterns without the need for stir bars. For SPME headspace and immersion sampling the special design (pat. pending) ensures optimal performance.

The PAL Heatex Stirrer offers:

- Rapid equilibration through effective stirring for headspace and immersion SPME sampling while ensuring the integrity of the fiber
- Efficient dissolution of solids, temperature controlled
- Thorough liquid/liquid extraction
- Stirring/heating for derivatization reactions
- No stir bar required, constant stirring also with samples containing solids
- No cross contamination
- Precise control of the equilibration temperature 40-150 °C
- Software controlled, temperature and stirring speed are logged

PAL SPME Arrow Ordering Information

System requirements


- PAL RTC or RSI with firmware 2.3 or higher
- PAL SPME Arrow Tool
- PAL Heatex Stirrer Module
- Adaptation of GC-injector (see page 11)
- A PAL SPME Arrow Conditioning Module is highly recommended, the PAL SPME Fiber Conditioning Module cannot be used with SPME Arrow.
- For overlapped extraction, the Agitator Module is required

The PAL SPME Arrows are available in order quantities of one, three or five SPME Arrows per box. For method development, a set of each fiber type (set of five) is available.

Diameter	Phase Thickness	Color Code	Set of 1 SPME Arrow Description PNo.	Set of 3 SPME Arrows Description PNo.	Set of 5 SPME Arrows Description PNo.		
PDMS SPME Arrow (Polydimethylsiloxane)							
1.1 mm	100 µm	Red	ARR11-P-100/20-P1	ARR11-P-100/20-P3	ARR11-P-100/20-P5		
Acrylate SPME Arrow (Polyacrylate)							
1.1 mm	100 µm	Grey	ARR11-A-100/20-P1	ARR11-A-100/20-P3	ARR11-A-100/20-P5		
Carbon WR SPME Arrow / PDMS (Carbon Wide Range / PDMS)							
1.1 mm	120 µm	Light Blue	ARR11-C-WR-120/20-P1	ARR11-C-WR-120/20-P3	ARR11-C-WR-120/20-P5		
DVB SPME Arrow / PDMS (Divinylbenzene / PDMS)							
1.1 mm	120 µm	Violet	ARR11-DVB-120/20-P1	ARR11-DVB-120/20-P3	ARR11-DVB-120/20-P5		
PDMS SPME Arrow (Polydimethylsiloxane)							
1.5 mm	250 µm	Black	ARR15-P-250/20-P1	ARR15-P-250/20-P3	ARR15-P-250/20-P5		
SPME Arrow Collection – Development Kit							
One of each t	ARR1115-SEL5-S1						

(Detailed user information can be found in the SPME Arrow Instruction Leaflet.)



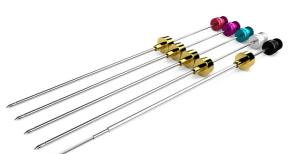

Ordering information for required parts



Ordering information for optional modules



Ordering information for starter kits and suitable liners

Starter Kits	PAL3-ARR-Start-GC2010	Starter Kit SPME Arrow less injector of Shimao Arrow Liner,1 SPME A
	PAL3-ARR-Start-GC6890	Starter Kit SPME Arrow injector of Agilent GC the Agilent GC 6890,
	PAL3-ARR-Start-GC7890	Starter Kit SPME Arrow injector of Agilent GC the Agilent GC 7890,
	PAL3-ARR-Start-Tr1300	Starter Kit SPME Arrow split/splitless injector o ers for SSL Injector of T
	PAL3-ARR-Start-TrUltra	Starter Kit SPME Arrow splitless injector of Th SSL Injector of Therm
Suitable Liners	ARR-Liner-CondModule	Liner for SPME Arrow
	ARR-Liner-GC2010	SPME Arrow Liner for
	ARR-Liner-GC6890	SPME Arrow Liner for
	ARR-Liner-GC7890	SPME Arrow Liner for
	ARR-Liner-Trace1300	SPME Arrow Liner for
	ARR-Liner-TraceUltra	SPME Arrow Liner for

tarter Kit containing SPME Arrow Tool SPME Arrow Holder SPME Arrow Selection of five different SPME Arrows SPME Arrow Instruction Leaflet Adapter Kit for GC injector (details see below)
or mixing and heating in sample prep and SPME Temperature range 30 - 150°C Mixing speed up to 1600 rpm (200 cycloidal loops) Optimized for 20 mL vials (for 10 mL vials a special adapter is required)
or the conditioning of SPME Arrow as well as SPME fibers rior to sample enrichment Position for automated conditioning Position for manual pre-conditioning Automated purge gas valve Manual gas valve for pre-conditioning Replacement liners for easy maintenance
he Agitator Module provides 6 positions for 20 mL vials for neubation and agitation of samples. Temperature range 40 – 200 °C Agitation speed 250 – 750 rpm Optional adapters for 2 mL or 10 mL vials
ow for Shimadzu GC 2010 consisting of: 1 Adaption Kit for the split/split- adzu GC 2010 (ARR-SSL-Inj-GC2010), 1 Liner Nut, 1 Screw Cap, 2 SPME Arrow Tool Kit (PAL3-SPME-Arrow-Kit)
ow for Agilent GC6890 consisting of: 1 Adaption Kit for the split/splitless C 6890 (ARR-SSL-Inj-GC6890), 2 SPME Arrow Liners for SSL Injector of 1 SPME Arrow Tool Kit (PAL3-SPME-Arrow-Kit)
ow for Agilent GC 7890 consisting of: 1 Adaption Kit for the split/splitless C 7890 (ARR-SSL-Inj-GC7890), 2 SPME Arrow Liners for SSL Injector of 1 SPME Arrow Tool Kit (PAL3-SPME-Arrow-Kit)
w for Thermo GC Trace 1300/1310 consisting of: 1 Adaptation Kit for the of Thermo GC Trace1300/1310 (ARR-SSL-Inj-Trace1300), 2 SPME Arrow Lin-Thermo GC Trace1300/1310, 1 SPME Arrow Tool Kit (PAL3-SPME-Arrow-Kit)
ow for Thermo GC Trace Ultra consisting of: 1 Adaptation Kit for the split/ hermo GC TraceUltra (ARR-SSL-Inj-TraceUltra), 2 SPME Arrow Liners for no GC TraceUltra, 1 SPME Arrow Tool Kit (PAL3-SPME-Arrow-Kit)
v Conditioning Module, package containing 3 pcs
r SSL Injector of SHIMADZU GC 2010, package containing 3 pcs
r SSL Injector of AGILENT GC 6890, package containing 3 pcs

r SSL Injector of AGILENT GC 7890, package containing 3 pcs

r SSL Injector of Thermo GC Trace1300, package containing 3 pcs

r SSL Injector of Thermo GC TraceUltra, package containing 3 pcs





Distributed by:

For more information on the PAL System visit:

www.palsystem.com



PAL is a registered trademark of CTC Analytics AG, Switzerland